


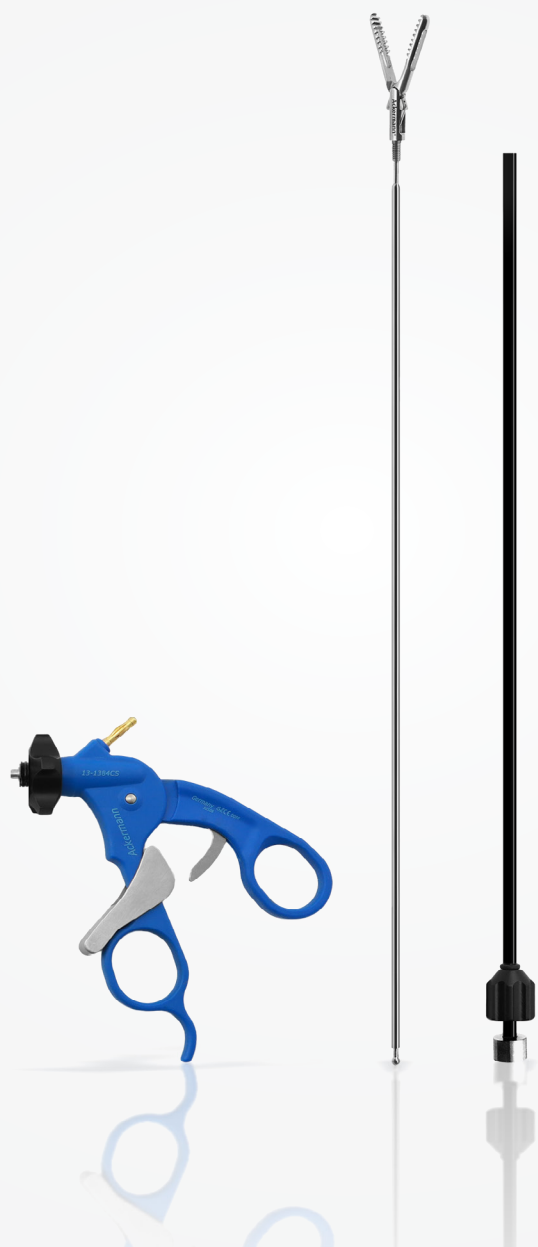


System CS™

Reusable Monopolar Instruments

System CS™ Line

Key Features

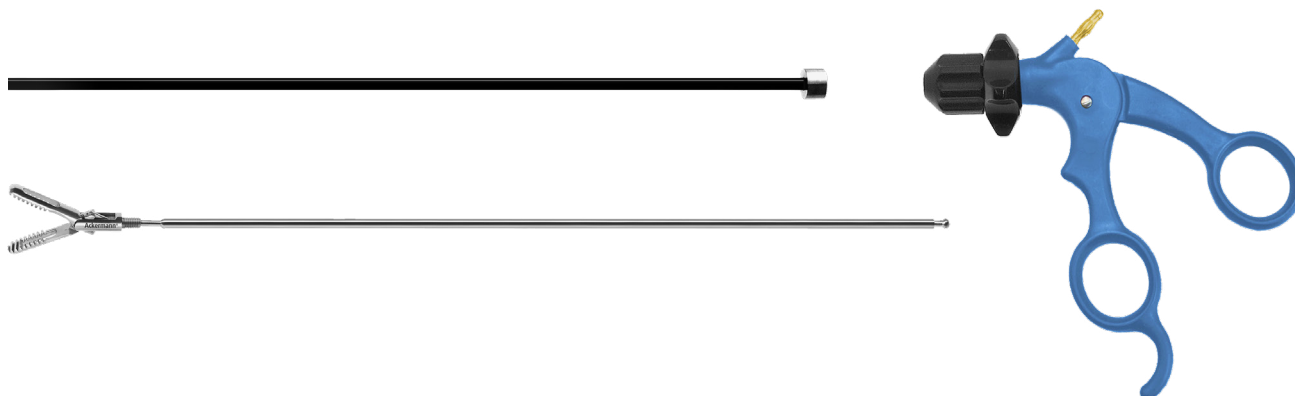


- Extremely easy assembly and disassembly with **screwing-system**
- **Premium performance** at unbeatable value for money without having to compromise on quality
- Comprehensive array of **instrument patterns**
- Wide selection of **handle options**

How the System works

Disassembly of System CS™ instruments

- 1 Open the handle completely and unscrew the nut counter clockwise.
- 2 Detach the insert/sheath from the handle by removing the ball coupling of the insert from the socket of the handle.
- 3 Unscrew the insert from the sheath and separate them fully.

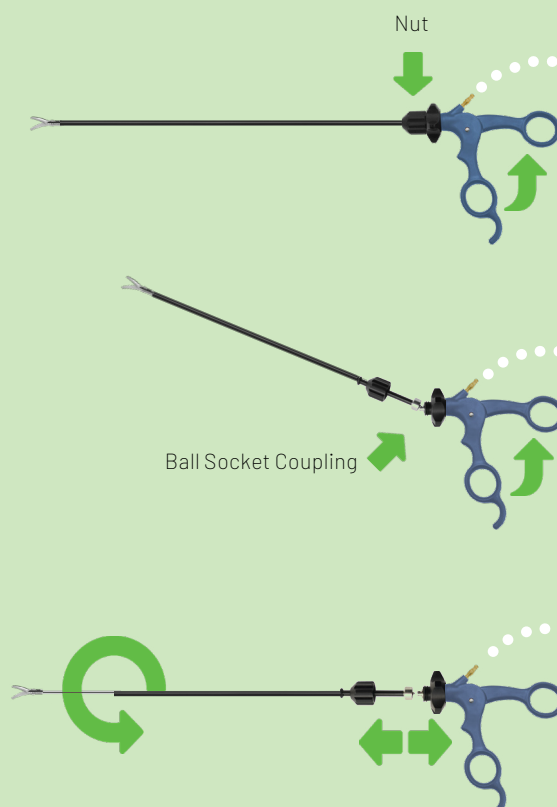



Wide Selection of Handle Options


Reference Identification

Select the desired handle type by numerical reference number:

- **#13-1367** Handle without ratchet
- **#13-1368** Handle with traditional ratchet
- **#13-1384** Handle with universal ratchet
- **#13-1369** Handle with french ratchet
- **#13-1358AXDI** Axial handle with ratchet

Select the desired color-coding for each System CS™ ring handle by adding the alphabetic characters:

- **CS**
- **FRIGELB**
- **FRIGRUEN**
- **FRIROT**

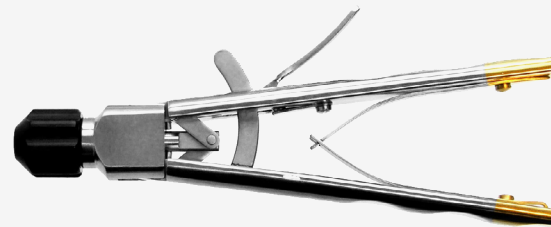

System CS™ Topsellers

Ordering Information and Recommendations

1x Metzenbaum scissors, curved, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle w/o ratchet

13-1461CSI

13-1356CS

13-1367CS



1x Hook scissors, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle w/o ratchet

13-1431CSI

13-1356CS

13-1367CS



1x Maryland dissector, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle w/o ratchet

13-1318CSI

13-1356CS

13-1367CS



1x Endoscopic clinch, atraumatic, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle with universal ratchet

13-1323CSI

13-1356CS

13-1368CS



1x Standard grasper, cross-serrated, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle with universal ratchet

13-1431CSI

13-1356CS

13-1384CS



1x DeBakey grasper, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle with universal ratchet

13-1570CSI

13-1356CS

13-1384CS



1x Johan grasper, slightly flat jaws, Ø 5 mm, 330 mm (*insert*)

1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm

1x System CS™ handle with universal ratchet

13-1312CSI

13-1356CS

13-1384CS




- 1x Babcock grasper, curved, ceramic, Ø 5 mm, 330 mm (*insert*)
- 1x Sheath for System CS™ inserts, Ø 5 mm, 330 mm
- 1x System CS™ handle with universal ratchet

13-1541CSI
13-1356CS
13-1384CS




- 1x Claw forceps, traumatic jaws, Ø 10 mm, 330 mm (*insert*)
- 1x Sheath for System CS™ inserts, Ø 10 mm, 330 mm
- 1x System CS™ handle with universal ratchet

13-1605CSI
13-1357CS
13-1384CS

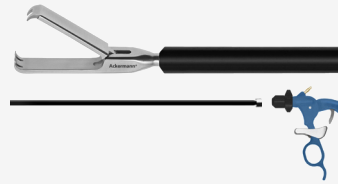

- 1x Biopsy spoon forceps, Ø 10 mm, 330 mm (*insert*)
- 1x Sheath for System CS™ inserts, Ø 10 mm, 330 mm
- 1x System CS™ handle w/o ratchet

13-1660CSI
13-1357CS
13-1367CS

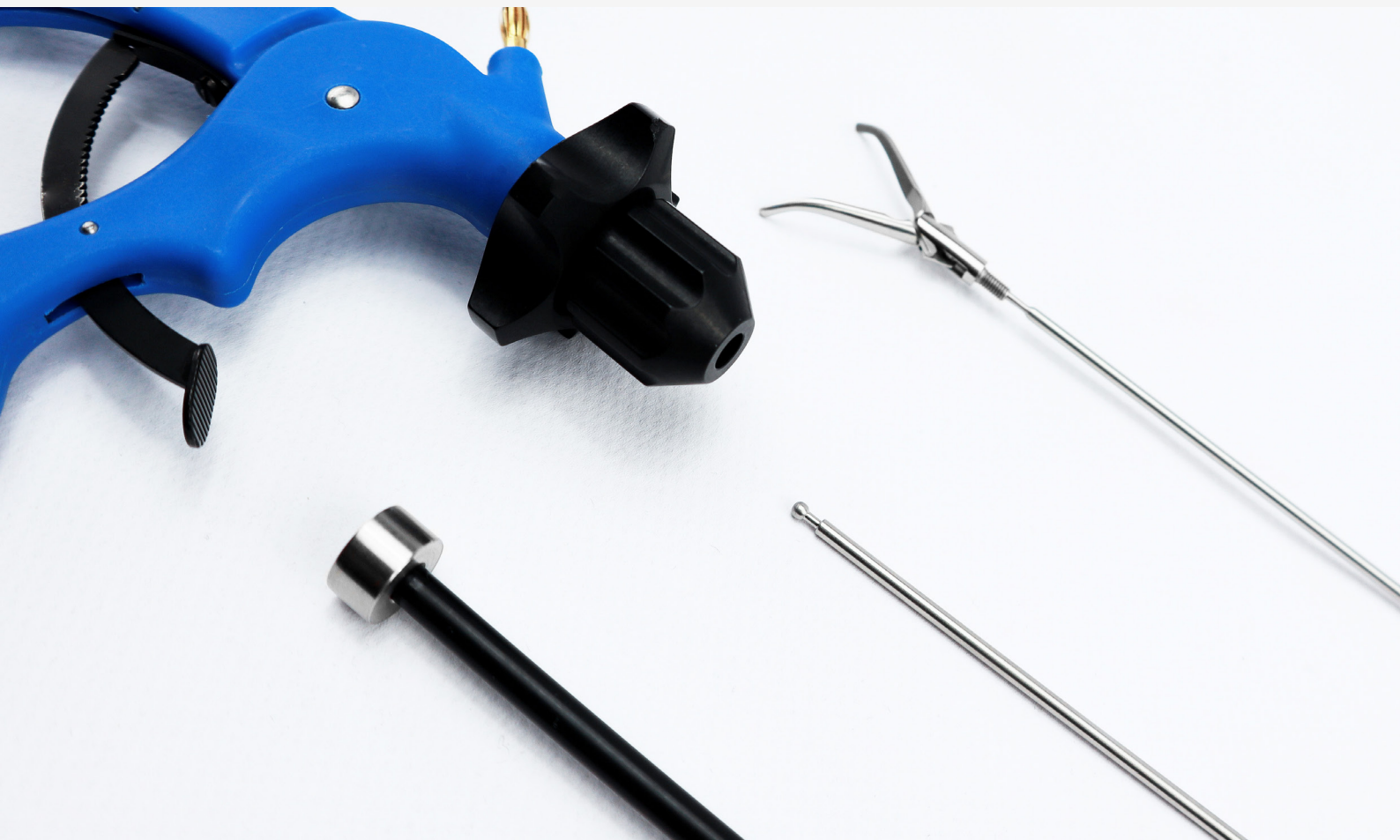



www.ackermannsurgical.com



Ackermann Instrumente GmbH
Eisenbahnstrasse 65-67
78604 Rietheim-Weilheim, Germany

Revision 04/22